



# U.L.E.M.A. in 8 steps (in a nutshell):


1. **Extract the demo data on your hard-disk.** The available data (*DemoDataUL* / *C3D*) are obtained from a typically developing child (left arm). The session folder contains all the C3D calibration files (each composed by the prefix *cal* and the name of the calibrated anatomical landmark) and 2 different tasks: HTH (Hand-to-Head) and RF (Reach Forward). There are 3 files per task (3 dynamic recordings), with 2 motion cycles (repetitions) per file. Each motion cycle is defined by 2 events (start and endpoint of the movement) of the type "General Event" (Context: General).
2. **Download and run the U.L.E.M.A. application.** After starting MATLAB, set the U.L.E.M.A. directory as the current working directory and type *mainGUI* in the Command window. The main window of the GUI should appear in a few seconds, after automatically including all the necessary paths in MATLAB.

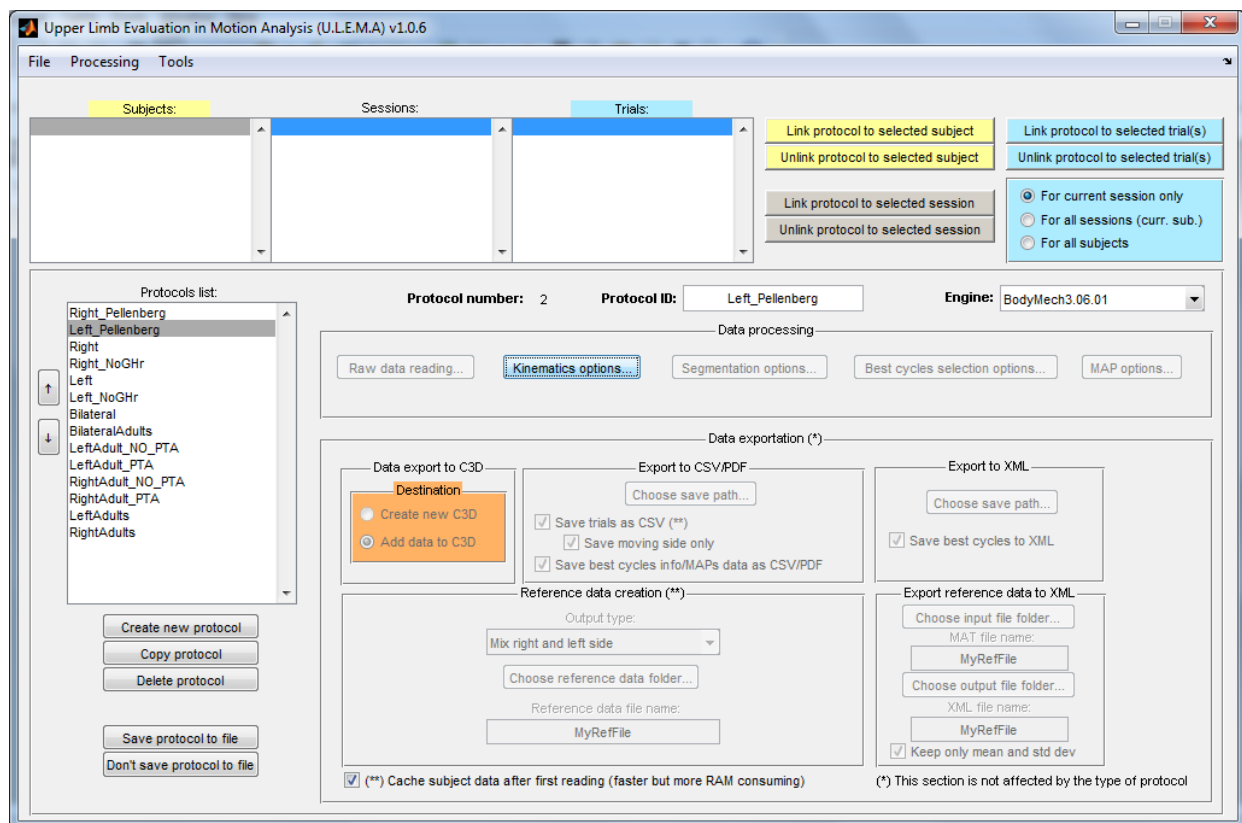

3. **Load the folder structure in the application.** Click on *File* → *Load* → *C3D tree (single subject)* and select the demo data folder *Subject1*. The structure should be loaded in the GUI.

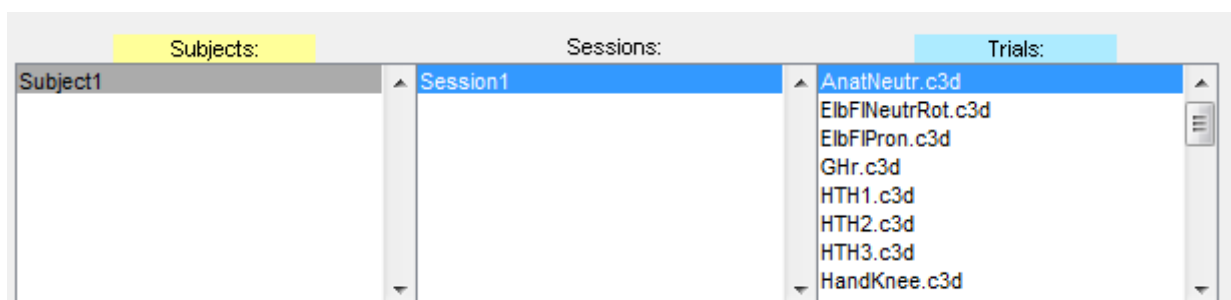


4. **Create a protocol with the correct processing settings.** A *protocol* in this context refers to a set of processing options that can manually be selected. The user can create a custom protocol by clicking on *Create new protocol*. A new item should be added to the list of available protocols. You can also change the name of the protocol (e.g. *DemoDataLeft*) by editing the field *Protocol ID*. After that, you should click on *Save protocol to file*. Modifications will be saved to file *ProtDB.mat* into the folder *Prot*. By clicking on the menu *Processing* → *Options* you can enable/disable single parts of the whole processing/exportation pipeline. If the section is checked in the menu, the corresponding button/panel will be activated in the main window. This will allow to modify the related processing options. Enable these sections and set the following parameters:
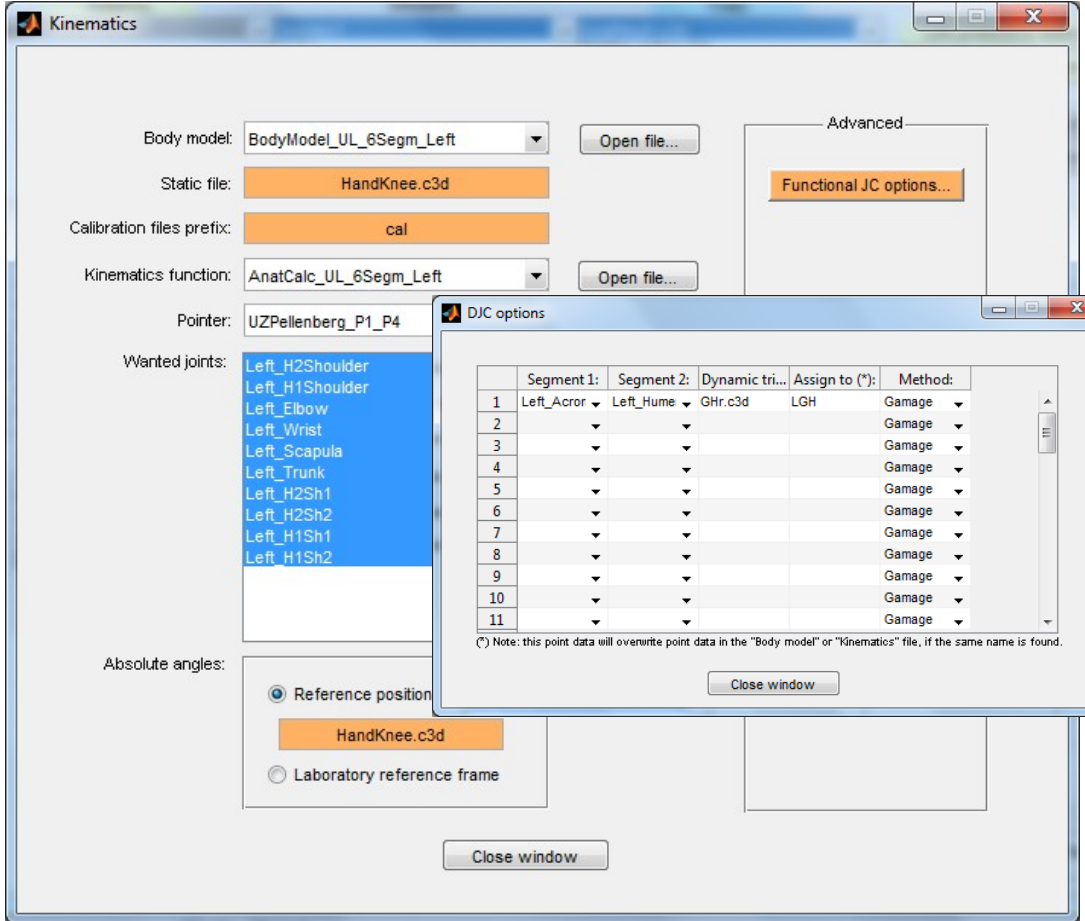
| Section       | Options  |
|---------------|--|
| Raw data read | There are no options (future releases), enable this section.   |
| Kinematics    |  <p><b>Body model (*)</b>: the BodyMech file defining the body model (in the folder <i>BodyModels</i>).</p> <p><b>Static file</b>: the static recording that will be used to define the reference position of the technical coordinate frames.</p> <p><b>Calibration files prefix</b>: string that must be present in the name of each calibration file (e.g. <i>cal</i>); the remaining part of the file name indicates the calibrated anatomical landmark.</p> <p><b>Kinematics function (*)</b>: the BodyMech file for calculating the joint angles (in the folder <i>AnatCalcs</i>). <b>Pointer (**)</b>: the BodyMech file for calculating the location of the calibrated anatomical landmarks (in the folder <i>Pointers</i>).</p> <p><b>Wanted joints (*)</b>: list of joints that are present in the selected body model file. Only the selected ones will be processed.</p> <p><b>Reference position (*)</b>: in a Body model file, some joints can have the proximal segment set at 0 (no explicit segment). For these joint, you can set the proximal reference frame as: <i>Laboratory reference frame</i>; a non-moving segment having the same 3D attitude of the distal segment in a specific file (<i>HandKnee.c3d</i> in the example).</p> |

**Functional Dynamic Joint Center:** every row represents the settings for a single joint centre to be computed. *Segment 1* and *2* are the proximal and distal segments respectively; *dynamic trial* refers to the trial containing the functional motion (i.e shoulder circumbduction); *assign to* indicates the name of the joint centre; *method* defines which method has to be used to compute the joint centre.

Segmentation

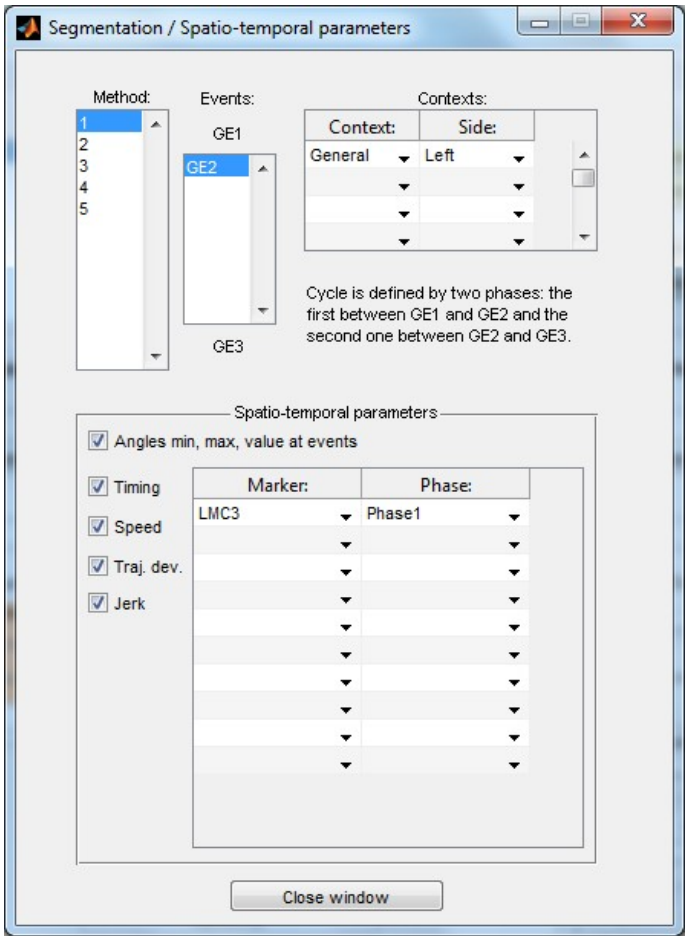

**Method:** method used for the segmentation; the full list can be found in the file *EventsConfig.txt* in the folder *Seg* (customisable). A brief description of the method is shown as well.

**Events:** list of events defining a motion cycle; the first and the last are used to split the entire data into separate motion cycles. In the example, *GE* stands for *General Event*.

**Contexts:** associates an event context (see Vicon Nexus) to a specific body side. This association can be done for multiple contexts.

**Spatiotemporal parameters:** each row of the table indicates a combination of marker and temporal phase in which the selected spatiotemporal parameters should be calculated. For every motion cycle, *Phase1* represents the time window between the first and second event (defined in Vicon Nexus), *phase2* the time between the second and third event, and so on. *EntireCycle* represents the entire motion cycle, even if multiple events are marked in the c3d file.

**Best cycles selection**

Task table:

| Task prefix: | Context:  | Phase:   | Edit angles:  |
|--------------|-----------|----------|---------------|
| HTH          | General ▼ | Phase1 ▼ | Click here... |
| HTM          | General ▼ | Phase1 ▼ | Click here... |
| HTS          | General ▼ | Phase1 ▼ | Click here... |
| RF           | General ▼ | Phase1 ▼ | Click here... |
| RGS          | General ▼ | Phase1 ▼ | Click here... |
| RGV          | General ▼ | Phase1 ▼ | Click here... |
| RS           | General ▼ | Phase1 ▼ | Click here... |
| RU           | General ▼ | Phase1 ▼ | Click here... |
|              | ▼         |          | Click here... |
|              | ▼         |          | Click here... |
|              | ▼         |          | Click here... |
|              | ▼         |          | Click here... |
|              | ▼         |          | Click here... |
|              | ▼         |          | Click here... |
|              | ▼         |          | Click here... |
|              | ▼         |          | Click here... |
|              | ▼         |          | Click here... |
|              | ▼         |          | Click here... |
|              | ▼         |          | Click here... |
|              | ▼         |          | Click here... |
|              | ▼         |          | Click here... |
|              | ▼         |          | Click here... |
|              | ▼         |          | Click here... |
|              | ▼         |          | Click here... |
|              | ▼         |          | Click here... |

Fast fill table ...

Angles list:

- LElevationplaneH1
- LShoulderElevationH1
- LShoulderRotationH1
- LEIbowFIExt
- LEIbowProSupination
- LWristFIExt
- LWristDeviation
- LScapProRetraction
- LScapRotation
- LScapTilting
- LTrunkFIExt
- LTrunkLateralFI
- LTrunkAxialRotation

Best cycles number:

Close window

**Best cycles selection**

**Task table:**

| Task prefix: | Context:  | Phase:   | Edit angles:  |
|--------------|-----------|----------|---------------|
| HTH          | General ▼ | Phase1 ▼ | Click here... |
| HTM          | General ▼ | Phase1 ▼ | Click here... |
| HTS          | General ▼ | Phase1 ▼ | Click here... |
| RF           | General ▼ | Phase1 ▼ | Click here... |
| RGS          | General ▼ | Phase1 ▼ | Click here... |
| RGV          | General ▼ | Phase1 ▼ | Click here... |
| RS           | General ▼ | Phase1 ▼ | Click here... |
| RU           | General ▼ | Phase1 ▼ | Click here... |
|              | ▼         |          | Click here... |
|              | ▼         |          | Click here... |
|              | ▼         |          | Click here... |
|              | ▼         |          | Click here... |
|              | ▼         |          | Click here... |
|              | ▼         |          | Click here... |
|              | ▼         |          | Click here... |
|              | ▼         |          | Click here... |
|              | ▼         |          | Click here... |
|              | ▼         |          | Click here... |
|              | ▼         |          | Click here... |
|              | ▼         |          | Click here... |
|              | ▼         |          | Click here... |
|              | ▼         |          | Click here... |
|              | ▼         |          | Click here... |




**Angles list:**

- LElevationplaneH1
- LShoulderElevationH1
- LShoulderRotationH1
- LEIbowFIExt
- LEIbowProSupination
- LWristFIExt
- LWristDeviation
- LScapProRetraction
- LScapRotation
- LScapTilting
- LTrunkFIExt
- LTrunkLateralFI
- LTrunkAxialRotation

Fast fill table ...


Best cycles number:

Close window

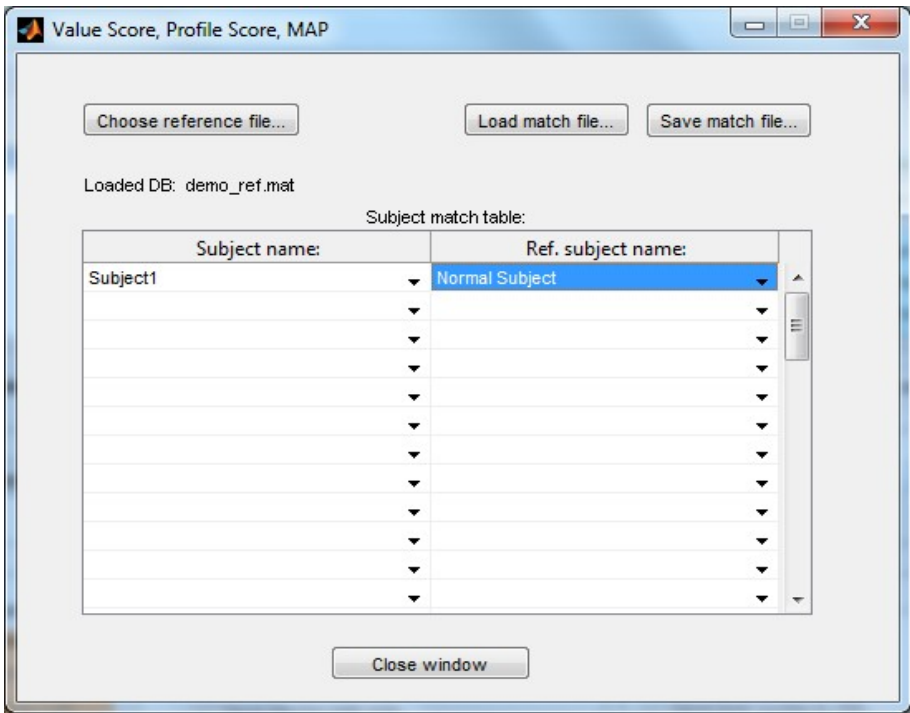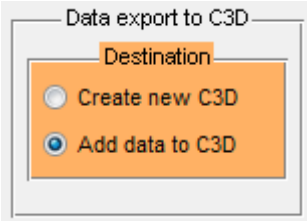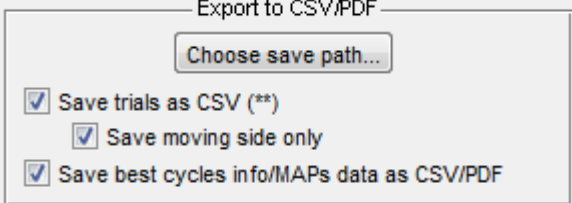
|                   |  |
|-------------------|--|
| MAP               |  <p><u>Choose reference file</u>: select the reference file to be used for the MAP calculation (reference database of typically developing children).</p> <p><u>Subject match table</u>: this table was once created to individually match a measured subject (patient) with a control subject. Currently, the reference database file created by U.L.E.M.A. contains only one subject (<i>Normal Subject</i>), which is the average based on 20 typically developing children. Hence, the only possible association is between each individual subject (patient) and this average normal subject.</p>   |
| Export to C3D     |  <p><u>Destination</u>: choose if you wish to create a new c3d file containing a copy of the original data of the loaded c3d files and all calculated data (calibrate anatomical landmarks, angles, ...): <i>Create new C3D</i>; or choose to add newly calculated data to the loaded ones (<i>Add data to C3D</i>). In the first case, a new subfolder will be created into the session folder containing the new c3d files.</p> <p><i>For this exportation option, it is mandatory to click on the menu Processing → Options → Set recovery path and select the folder containing the original data structure. In this example, that would be the folder containing the subfolder Subject1. The software is now able to automatically retrieve the original C3D files, copy and/or append data to them, and save them again.</i></p> |
| Export to CSV/PDF |  <p><u>Choose save path</u>: choose the folder where the CSV files (containing waveform data, best cycles ranking data...) and the PDF files (containing the MAP graphs) will be saved.</p> <p><u>Save moving side only</u>: if checked, each exported motion cycles will only contain data for angles belonging to the same side as the cycle. For instance, for a right side cycle, only angles whose name begins with R will be exported.</p>   |

(\*) Read BodyMech documentation for more details. (\*\*) The pointer is a custom made tool (with concurrent code) for every centre.

5. **Link the protocol to the selected trials.** This way, the software knows how to process each of the C3D files. Each file (dynamic trial) can have only one linked protocol. However, multiple different protocols can be used within one session, e.g. different protocols for different tasks. Select *DemoDataLeft* protocol and the trials RF1.c3d, RF2.c3d, RF3.c3d, HTH1.c3d, HTH2.c3d, HTH3.c3d, then click on *Link protocol to selected trial(s)*. The protocol name should be appended to the right of the files name.

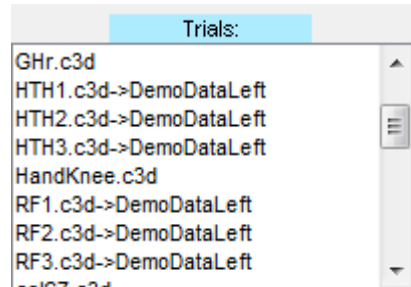

6. **Process data (finally!).** Before the processing can take place, click on *Processing* → *Options* → *Set proc. data path*. and select a folder that will contain the MAT files created during data processing. You can now start the processing phase by clicking on *Processing* → *RUN*. For U.L.E.M.A. 1.0.6+, you must also select the recovery path (*Processing* → *Options* → *Set recovery path*): choose *DemoDataUL / C3D*.
7. **Check the new content of into C3D files.** After exporting new data to C3D, reconstructed anatomical landmarks and joint angles can be visualised with your preferred software ([Vicon Nexus](#) / [Polygon](#), [Mokka](#),...) as any other marker. Spatiotemporal parameters are also stored into the C3Ds. Below is an example using Vicon Nexus (reconstructed anatomical landmarks are depicted in grey, while the technical markers are depicted in yellow):

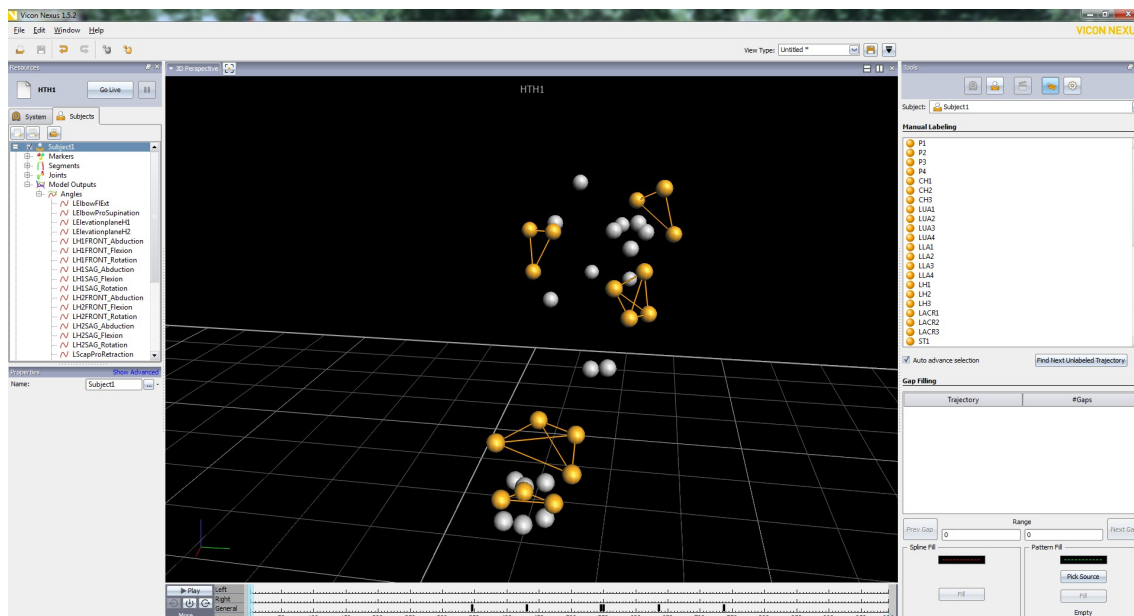

8. **Wait! I would like to change some processing settings and re-process some sections. Can I do that?.** Of course! You can load the MAT files that were saved during the processing phase (in the processing data path). By clicking on *File* → *Load* you can choose to load one single MAT file or many MAT files in batch (read full documentation). And you can re-process whatever you want (don't forget to enable the option *Processing* → *Options* → *Force re-processing*).

For any detail, please read the full doc PDF contained in every distribution of the application.